



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D. F-ABFT  
Chairman

Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

Joshua Smith  
Director

ATTESTATION

I, the undersigned Joshua Smith (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

- 1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-003314, in compliance with the agency's Maintenance Policy and Procedures.

Signed

Joshua Smith

Name of Position: Director/Records Custodian


State of Oklahoma

County of Oklahoma

Signed or attested before me on 06/01/2026 by

Joshua Smith.

(Signature of notarial officer)

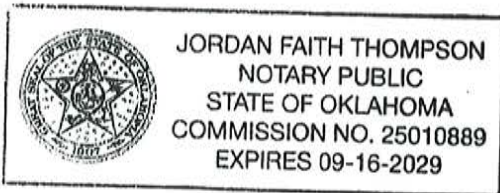

# Service Overview

Board of Tests for Alcohol and Drug Influence  
Enhancing public safety by ensuring the accuracy and scientific reliability  
of blood and breath testing



3600 N. Martin Luther King Bldg #1  
Oklahoma City, OK 73111  
405-425-2460

**BENCH CHECK DATE: 5/29/2026**

**SERIAL: 80-003314**

## **DOCUMENTATION WITHIN SERVICE PACKET:**

- BOT Certificate of Maintenance, Calibration, and Operation
- Bench Check
- BOT Dry Gas Certificate of Analysis
- NHTSA Supplemental CPL Letter
- Intermountain Specialty Gases Certificate of Analysis
- Mock Subject Test
- BOT Certificate of Analysis on Certified Solutions (4)
- Guth Laboratories Certificate of Analysis (4)



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

**CERTIFICATE of MAINTENANCE, CALIBRATION, and OPERATION**

This certifies the Intoxilyzer 8000 referenced below, manufactured by CMI, Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, is maintained by the Board of Tests for Alcohol and Drug Influence in accordance with 40 O.A.C. 30-1-3. The instrument reported below is hereby deemed to be in compliance with Board rules via linearity accuracy examination and operational/functional procedure evaluation also referred to as a "Bench Check". The Intoxilyzer 8000 appears on the National Highway Traffic Safety Administration's Conforming Products List of Evidential Breath Alcohol Measurement Devices (82 FR 50940 *et seq*) as an approved device. In accordance with 40 O.A.C. 25-1-2, this instrument is approved and certified for use in the State of Oklahoma as an automated evidential breath-alcohol measurement device by the Board of Tests pursuant to 47 O.S. § 759.

All analyses conducted by this instrument are considered precise, valid, and in accordance with 40 O.A.C. 30-1-3(b) Operating Procedure(s). A printed breath test operator affidavit containing the seal of the Board by this instrument confirms a calibration check successfully verified the instrument's calibration at the time of test. The printed breath test affidavit also serves as proof that safeguards in place to protect the breath samples, including but not limited to Ambient Air Interference, Radio Frequency Interference, Residual Mouth Alcohol, Deficient Samples, etc. were not triggered and the test is a true and accurate scientific measurement of the person's alcohol concentration.

- Maintenance was performed on the Intoxilyzer 8000, Serial No. 80-003314 on 5/29/2026 in accordance with 40 O.A.C. 30-1-3 and internal BOT maintenance policy and procedures.
- A dry gas canister appearing on the National Highway Traffic Safety Administration's Conforming Products List of Calibrating Units for Breath Alcohol Testers (77 FR 64588 *et seq*) and approved in accordance with 40 O.A.C. 25-1-3 and BOT maintenance policy was installed.
  - Manufacturer: Intermountain Unit: 208 ppm / 105 L
  - Cylinder Lot No. - Canister No.: 26-4437-195
  - BOT analysis average test result: 0.081
  - Cylinder expiration date: 1/29/2028
- In accordance with 40 O.A.C. 25-1-2.1, the approved reference standards and equipment were used:

| Simulator Model          |             |             |             |
|--------------------------|-------------|-------------|-------------|
| Guth 12V500              | Guth 12V500 | Guth 12V500 | Guth 12V500 |
| Serial #                 |             |             |             |
| MP9201                   | MP9202      | MP9203      | MP9204      |
| Concentration            |             |             |             |
| .02                      | .04         | .10         | .20         |
| Lot #                    |             |             |             |
| 25310                    | 25400       | 25282       | 25360       |
| Manufactured Date        |             |             |             |
| 9/3/2025                 | 12/3/2025   | 8/12/2025   | 10/21/2025  |
| Expiration Date          |             |             |             |
| 9/3/2027                 | 12/3/2027   | 8/12/2027   | 10/21/2027  |
| Solution Commission Date |             |             |             |
| 5/28/2026                | 5/28/2026   | 5/28/2026   | 5/28/2026   |

Completed by BOT Personnel:

  
Signature \_\_\_\_\_  
Name Christopher Pape

Breath-Alcohol Testing Program Administrator  
Title \_\_\_\_\_

  
Signature \_\_\_\_\_  
Name Joshua Smith

State Director  
Title \_\_\_\_\_

## STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT

|                                       |  |  |
|---------------------------------------|--|--|
| BENCH CHECK DATE<br><b>05/29/2026</b> | TIME PROCEDURE BEGAN<br><b>09:52:0</b> | INSTRUMENT LOCATION<br><b>80-003314</b>      |
| INSTRUMENT TYPE<br><b>Intoxilyzer</b> | INSTRUMENT MODEL<br><b>I-8000</b>      | INSTRUMENT SERIAL NUMBER<br><b>80-003314</b> |

### TEST RESULTS

**Diagnostics**  
PASS

**Mouth Alcohol Test**  
PASS

**RFI Test**  
PASS

**Abort Test**  
PASS

**Improper Sample Test**  
PASS

**Interferent Detect Test**  
PASS

**No Sample Given Test**  
PASS

**Insufficient Sample Test**  
PASS

**Air Blank Contamination Test**  
PASS

**Range Exceeded Test**  
PASS

#### Duplicate Breath Sampling Results

| Value | Test 1 | Test 2 | Control Test |
|-------|--------|--------|--------------|
| .02   | 0.019  | 0.020  | 0.020        |
| .04   | 0.040  | 0.040  | 0.041        |
| .10   | 0.100  | 0.099  | 0.102        |
| .20   | 0.198  | 0.198  | 0.202        |

PASS

### CERTIFYING TECHNICIAN

I swear under penalty of perjury that in regards to the above listed instrument, I have complied with the Intoxilyzer Bench Checking Procedure approved by the State of Oklahoma.

|  |                               |
|--|-------------------------------|
| SIGNATURE<br> | DATE<br><b>5/29/2026</b>      |
| NAME<br><b>PAPE, CHRISTOPHER</b>   | PERMIT #<br><b>0000023020</b> |



U.S. Department  
of Transportation  
**National Highway  
Traffic Safety  
Administration**

1200 New Jersey Avenue SE.  
Washington, DC 20590

October 24, 2022

Scott Feely  
Intermountain Specialty Gases  
21913 Cobalt Ave.  
Caldwell, ID 83605

Dear Mr. Feely:

This letter confirms that Intermountain Specialty Gases submitted the following model dry gas calibration device filled a fixed-concentration of ethanol-in-nitrogen was tested:

- **Model: 105S (50 ppm-550ppm)**

This calibrating unit previously conformed to the Model Specifications and is listed on the Conforming Products List of Alcohol Calibration Units. The device is being reevaluated because Intermountain Specialty Gases relocated to a new filling facility. The cylinder also has a new DOT rating that allows flammable gas. This is simply a rating change and involved no change to the cylinder.

The ethanol-in-nitrogen dry gas standards were evaluated by the Volpe National Transportation Systems Center and found to conform to the model specifications for calibrating units (CUs) for breath alcohol testers. These tanks will be included in the next update of the Conforming Products List for CUs. A copy of the evaluation report is provided for your information.

NHTSA has also reviewed the Quality Assurance Plan (QAP) for the **Intermountain Specialty Gases 105S** ethanol-in-nitrogen gas standards dated February 3, 2022, as required by 49 CFR Part 40, "Procedures for Transportation Drug and Alcohol Testing Program" (59 FR 7340). We are pleased to inform you that this submission has been approved.

We will place this QAP on file at NHTSA, together with a copy of this letter for interested agencies. A copy will also be provided to the U.S. Department of Transportation Office of Drug and Alcohol Policy and Compliance for their records.

Please be advised that this office is to be notified immediately of any revisions or modifications to the ethanol-in-nitrogen standards or this QAP.

If you require further assistance or have any questions regarding NHTSA's quality assurance program in relation to 49 CFR Part 40, please do not hesitate to contact me at [NHTSA.Alcohol@dot.gov](mailto:NHTSA.Alcohol@dot.gov).

Sincerely,



DeReece Smither, Ph.D.  
Research Psychologist  
Behavioral Research Division

Enclosure

CC: Patrice Kelly, OST  
Bohdan Baczara, OST  
Edward Conde, Volpe



# STATE OF OKLAHOMA BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Dr. Jarrad Wagner, Ph.D., F-ABFT  
Chairman

Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

Joshua Smith  
Director

## Certificate of Analysis – Dry Gas

ALPHA INSTRUMENT  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-002591  
04/27/2026 11:26

| Test            | g/210L | Time  |
|-----------------|--------|-------|
| Air Blank       | 0.000  | 11:26 |
| Cal Check       | 0.081  | 11:27 |
| Air Blank       | 0.000  | 11:27 |
| Cal Check       | 0.080  | 11:27 |
| Air Blank       | 0.000  | 11:28 |
| Cal Check       | 0.081  | 11:28 |
| Air Blank       | 0.000  | 11:29 |
| Cal Check       | 0.081  | 11:29 |
| Air Blank       | 0.000  | 11:30 |
| Cal Check       | 0.081  | 11:30 |
| Air Blank       | 0.000  | 11:31 |
| Cal Check       | 0.081  | 11:31 |
| Air Blank       | 0.000  | 11:31 |
| Cal Check       | 0.081  | 11:32 |
| Air Blank       | 0.000  | 11:32 |
| Cal Check       | 0.081  | 11:33 |
| Air Blank       | 0.000  | 11:33 |
| Cal Check       | 0.081  | 11:33 |
| Air Blank       | 0.000  | 11:34 |
| Cal Check       | 0.081  | 11:34 |
| Air Blank       | 0.000  | 11:35 |
| Cal Check Stats |        |       |
| Average         | 0.0809 |       |
| Std Dev         | 0.0003 |       |
| Rel Std Dev (%) | 0.3908 |       |

| Known Breath-Alcohol Equivalent (g/210L) Ratio of Nitrogen and Ethanol Gas (208 ppm/105 L)<br>0.08 g/210L = 208 PPM + 2605*   |  |                                    |
|---|--|------------------------------------|
| Manufacturer: <b>Intermountain</b>  | Manufacturer Cert. No.: <b>C252001</b> | Cylinder Size: <b>105L</b>         |
| Lot No.: <b>26-4437</b>   | Cylinder No.: <b>195</b>               | Expiration Date: <b>01/29/2028</b> |
| The above 105 L cylinder has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference standard meets the requirements for use as a calibrating unit in accordance with 47 O.S. 759, 40 OAC 25-1-3, and appears on the current conforming products list of calibrating units for breath alcohol testers published by the U.S. Department of Transportation in the Federal Register [77 FR 64588 et seq]. * See Dubowski, K.M. and Essary N., <i>Journal of Anal. Toxicol.</i> (1996) 20(6): pp. 484-491. [g/210L = PPM +2605]. |  |                                    |
| Signature<br>   | Date:<br><b>04/27/2026</b>             |                                    |
| Name<br><b>Christopher Pape</b>   | Permit No.: <b>23020</b>               |                                    |



Idaho H.Q.  
21913 Cobalt Ave.  
Caldwell, ID 83605  
(208) 585-5829

Texas  
5041 Spencer Hwy. #707  
Pasadena, TX 77505  
www.isgases.com

Canada  
#4, 4830 - 78th St.  
Red Deer, AB T4P 3B2  
orders@isgases.com

Version 1.0

## CERTIFICATE OF ANALYSIS

~ for Calibration and Analytical Use ~

| Lot Number                | Part Number     | Manufacture Date | Expiration Date |
|---------------------------|-----------------|------------------|-----------------|
| 26-4437                   | 14-0208         | January 29, 2026 | January, 2028   |
|                           | Cylinder Volume | Pressure         | Valve           |
|                           | 105L            | 1,050 psi        | C-10            |
| Parent Cylinder ID Number | ND7020, ND7067  |                  |                 |

| Mixture Specifications |       |               |   |
|------------------------|-------|---------------|---|
| Components             | BrAC  | Certification | Uncertainty (relative) ( $\pm$ )                    |
| Ethanol                | 0.08% | 208.3 ppm     | 2% or 0.002 BrAC (g/210L),<br>whichever is greater. |
| Nitrogen               | —     | Balance       |   |

| Method of Preparation  | Preparation Accuracy ( $\pm 2\%$ ) |
|--|------------------------------------|
| This mix was prepared gravimetrically and is traceable to the NIST by certified weights (ID #CA10814) used to verify the calibration of the scale. The parent cylinder was analyzed against standards traceable to VSL* (the Dutch National Metrology Institute) on a FT-IR. |                                    |

| Primary Reference Material |                                   |                        |                |
|----------------------------|-----------------------------------|------------------------|----------------|
| Cylinder #                 | D245432                           | Certificate #          | C2552001       |
| Certified Date             | May 23, 2025                      | Expire Date            | April 30, 2028 |
| Certified Value            | 127.9 ppm Ethanol/ N <sub>2</sub> | Analytical Uncertainty | $\pm 1.0\%$    |

| Safety and Storage  |
|---|
| Only use product for intended use and with correct regulators/fittings. Handle product safely and do not puncture vessel. Store product in a secure area (with cap, when possible) and do not expose product to temperatures $< 45^{\circ}\text{F}$ ( $7^{\circ}\text{C}$ ) for an extended period of time. |

RCVD  
4/17/26  
CP

Results Authorized by: Aliya Burnside  
Aliya Burnside  
Title: Laboratory Manager  
Certificate Issue Date: January 29, 2026

\*VSL was deemed equivalent to the NIST in a memo issued July 1, 2022 (available upon request).

This certificate is to certify that the cylinders referenced have been tested and verified to meet the defined specifications. The tests were performed using gases and equipment that are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI). The basis of compliance stated is a comparison of the measurement parameters to the specified or required calibration process. The expanded uncertainties use a coverage factor of  $k=2$  to approximate the 95% confidence level of the measurement, unless otherwise noted. This certificate only applies to the components described and contained in this particular lot. This certificate shall not be reproduced except in full without written approval from Intermountain Specialty Gases (ISG). If not included, the uncertainty of calibrations/tests are available upon request. ISG is accredited to ISO/IEC 17025:2017 for testing and to ISO 17034:2016 for reference material and certified reference material production by PJLA, Accreditation No. 99538.

End of Document

# OFFICER'S IMPAIRED DRIVING AFFIDAVIT

9/19

Section 1

|                                       |  |                     |                               |                             |             |             |              |                    |                   |
|---------------------------------------|--|---------------------|-------------------------------|-----------------------------|-------------|-------------|--------------|--------------------|-------------------|
| ARREST DATE<br>05/29/2026             |  | ARREST TIME<br>1045 |                               | COUNTY #<br>55              |             | CITY #      |              | CITATION #         |                   |
| ARREST LOCATION<br>BOT                |  |                     |                               |                             | CITY        |             |              | COUNTY<br>OKLAHOMA |                   |
| SUBJECT NAME<br>NORTH, GAYLON LEVERNE |  |                     |                               | DATE OF BIRTH<br>08/09/1961 |             | SEX<br>M    |              | HEIGHT<br>510      | WEIGHT<br>165     |
| ADDRESS<br>3600 MARTIN LUTHER KING A  |  |                     |                               | CITY<br>OKLAHOMA CITY       |             |             | STATE<br>OK  |                    | ZIP CODE<br>73133 |
| DRIVER LICENSE #<br>T083454871        |  |                     | EXPIRATION DATE<br>06/30/2019 |                             | STATE<br>OK | CLASS<br>A  | RESTRICTIONS | ENDORSEMENTS<br>M  |                   |
| VEHICLE MAKE<br>MOCK                  |  | MODEL<br>MOCK       |                               | TAG #<br>MOCK               |             | STATE<br>OK | CMV<br>N     | PLACARD<br>N       |                   |

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST

Describe person's condition (odor, actions, etc.):

ROUTINE MAINTENANCE, INSTALLED NEW CYLINDER

**THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST**

**BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.**

Section 2

| EQUIPMENT<br>I-8000   |                         | SERIAL NUMBER<br>80-003314   |        | SPECIALIST<br>BOARD OF TESTS   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|---|-------------------------|--|--------|--|-------------------------------------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|------|--|-------|-----------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|-----------|------|-------|-----------|------|-------|
| STANDARD<br>Dry   |                         | TARGET VALUE<br>0.080  |        | PERMIT NUMBER<br>23020   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| MANUFACTURER<br>INTERMOUNTAIN   |                         | LOT #<br>26-4437   |        | DATE CYLINDER INSTALLED<br>05/29/2026  | DATE CYLINDER EXPIRES<br>01/29/2028 |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| OBSERVATION START<br>1050   | TEST DATE<br>05/29/2026 | <p>A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p>I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;"><i>Christopher Pape</i><br/>SIGNATURE OF OPERATOR</p> |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| RESULT TABLE  |                         |  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| <table style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="text-align: left;">Test</th> <th style="text-align: left;">g/210L</th> <th style="text-align: left;">Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>11:13</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:13</td></tr> <tr><td style="padding-left: 20px;">Breath Volume</td><td>1.968 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:14</td></tr> <tr><td>Wait</td><td></td><td>11:16</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:16</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:16</td></tr> <tr><td style="padding-left: 20px;">Breath Volume</td><td>1.871 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:17</td></tr> <tr><td>Cal Check</td><td>0.08</td><td>11:17</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:18</td></tr> </tbody> </table> |                         | Test   | g/210L | Time   | Air Blank                           | 0.00 | 11:13 | Subject Test | 0.00 | 11:13 | Breath Volume | 1.968 LITERS |  | Air Blank | 0.00 | 11:14 | Wait |  | 11:16 | Air Blank | 0.00 | 11:16 | Subject Test | 0.00 | 11:16 | Breath Volume | 1.871 LITERS |  | Air Blank | 0.00 | 11:17 | Cal Check | 0.08 | 11:17 | Air Blank | 0.00 | 11:18 |
| Test  | g/210L                  | Time   |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:13  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Subject Test  | 0.00                    | 11:13  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Breath Volume   | 1.968 LITERS            |  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:14  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Wait  |                         | 11:16  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:16  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Subject Test  | 0.00                    | 11:16  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Breath Volume   | 1.871 LITERS            |  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:17  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Cal Check   | 0.08                    | 11:17  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:18  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | NAME<br>PAPE, CHRISTOPHER  |        | <div style="border: 1px solid black; padding: 5px; text-align: center;"> <p>Test Result<br/>BrAC g/210L<br/><b>0.00</b></p> </div> |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | BADGE #<br>0000000484  |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | PERMIT #<br>0000023020   |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | AGENCY<br>BOARD OF TESTS   |        |  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |



Section 3

**What happens next?**

1. Upon receipt and review of the Officer's Impaired Driving Affidavit and any associated test results, if applicable, the State of Oklahoma may revoke your Oklahoma driving privileges.
2. If you hold a commercial driver license or were driving a commercial motor vehicle, the State of Oklahoma may disqualify your privilege to operate a commercial motor vehicle.
3. You will continue to have driving privileges, if otherwise eligible, unless or until an Order of Revocation and/or Disqualification is issued by the State of Oklahoma.
4. The State of Oklahoma will send you the Order of Revocation and/or Disqualification to the address you have on file with Service Oklahoma.  
To confirm or change your address on file, contact Service Oklahoma.
5. The Order and Notice of Revocation will inform you of your requirements to regain driving privileges in Oklahoma.

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: \_\_\_\_\_ Place (location when signed): \_\_\_\_\_ Signature of arresting officer: \_\_\_\_\_

|                           |                                |                          |                     |
|---------------------------|--------------------------------|--------------------------|---------------------|
| NAME<br>PAPE, CHRISTOPHER |                                | AGENCY<br>BOARD OF TESTS |                     |
| BADGE<br>0000000484       | AGENCY ADDRESS<br>INTER-AGENCY | ZIP CODE<br>00000        | PHONE<br>4054252460 |

**OTHER WITNESSES:**

Section 5

|   |      |       |         |       |
|---|------|-------|---------|-------|
| 1 | NAME | TITLE | ADDRESS | PHONE |
| 2 | NAME | TITLE | ADDRESS | PHONE |

# OFFICER'S IMPAIRED DRIVING AFFIDAVIT

10/19

Section 1

|                                       |  |                     |                               |                             |             |             |                    |                   |               |
|---------------------------------------|--|---------------------|-------------------------------|-----------------------------|-------------|-------------|--------------------|-------------------|---------------|
| ARREST DATE<br>05/29/2026             |  | ARREST TIME<br>1045 |                               | COUNTY #<br>55              |             | CITY #      |                    | CITATION #        |               |
| ARREST LOCATION<br>BOT                |  |                     |                               | CITY                        |             |             | COUNTY<br>OKLAHOMA |                   |               |
| SUBJECT NAME<br>NORTH, GAYLON LEVERNE |  |                     |                               | DATE OF BIRTH<br>08/09/1961 |             | SEX<br>M    |                    | HEIGHT<br>510     | WEIGHT<br>165 |
| ADDRESS<br>3600 MARTIN LUTHER KING A  |  |                     | CITY<br>OKLAHOMA CITY         |                             |             | STATE<br>OK |                    | ZIP CODE<br>73133 |               |
| DRIVER LICENSE #<br>T083454871        |  |                     | EXPIRATION DATE<br>06/30/2019 |                             | STATE<br>OK | CLASS<br>A  | RESTRICTIONS       | ENDORSEMENTS<br>M |               |
| VEHICLE MAKE<br>MOCK                  |  | MODEL<br>MOCK       |                               | TAG #<br>MOCK               |             | STATE<br>OK | CMV<br>N           | PLACARD<br>N      |               |

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST


Describe person's condition (odor, actions, etc.):

ROUTINE MAINTENANCE, INSTALLED NEW CYLINDER

**THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST**

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

| EQUIPMENT<br>I-8000   |                         | SERIAL NUMBER<br>80-003314   |        | SPECIALIST<br>BOARD OF TESTS  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|---|-------------------------|--|--------|---|-------------------------------------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|------|--|-------|-----------|------|-------|--------------|------|-------|---------------|--------------|--|-----------|------|-------|-----------|------|-------|-----------|------|-------|
| STANDARD<br>Dry   |                         | TARGET VALUE<br>0.080  |        | PERMIT NUMBER<br>23020  |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| MANUFACTURER<br>INTERMOUNTAIN   |                         | LOT #<br>26-4437   |        | DATE CYLINDER INSTALLED<br>05/29/2026   | DATE CYLINDER EXPIRES<br>01/29/2028 |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| OBSERVATION START<br>1050   | TEST DATE<br>05/29/2026 | <p>A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.</p> <p>I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.</p> <p style="text-align: center;"><i>Christopher Pape</i><br/>SIGNATURE OF OPERATOR</p> |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| RESULT TABLE  |                         |  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| <table style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="text-align: left;">Test</th> <th style="text-align: left;">g/210L</th> <th style="text-align: left;">Time</th> </tr> </thead> <tbody> <tr><td>Air Blank</td><td>0.00</td><td>11:13</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:13</td></tr> <tr><td style="padding-left: 20px;">Breath Volume</td><td>1.968 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:14</td></tr> <tr><td>Wait</td><td></td><td>11:16</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:16</td></tr> <tr><td>Subject Test</td><td>0.00</td><td>11:16</td></tr> <tr><td style="padding-left: 20px;">Breath Volume</td><td>1.871 LITERS</td><td></td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:17</td></tr> <tr><td>Cal Check</td><td>0.08</td><td>11:17</td></tr> <tr><td>Air Blank</td><td>0.00</td><td>11:18</td></tr> </tbody> </table> |                         | Test   | g/210L | Time  | Air Blank                           | 0.00 | 11:13 | Subject Test | 0.00 | 11:13 | Breath Volume | 1.968 LITERS |  | Air Blank | 0.00 | 11:14 | Wait |  | 11:16 | Air Blank | 0.00 | 11:16 | Subject Test | 0.00 | 11:16 | Breath Volume | 1.871 LITERS |  | Air Blank | 0.00 | 11:17 | Cal Check | 0.08 | 11:17 | Air Blank | 0.00 | 11:18 |
| Test  | g/210L                  | Time   |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:13  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Subject Test  | 0.00                    | 11:13  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Breath Volume   | 1.968 LITERS            |  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:14  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Wait  |                         | 11:16  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:16  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Subject Test  | 0.00                    | 11:16  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Breath Volume   | 1.871 LITERS            |  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:17  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Cal Check   | 0.08                    | 11:17  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
| Air Blank   | 0.00                    | 11:18  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | NAME<br>PAPE, CHRISTOPHER  |        | <br><div style="border: 1px solid black; padding: 5px; text-align: center;">                     Test Result<br/>BrAC g/210L<br/><b>0.00</b> </div> |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | BADGE #<br>0000000484  |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | PERMIT #<br>0000023020   |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |
|   |                         | AGENCY<br>BOARD OF TESTS   |        |   |                                     |      |       |              |      |       |               |              |  |           |      |       |      |  |       |           |      |       |              |      |       |               |              |  |           |      |       |           |      |       |           |      |       |

Section 3

**What happens next?**

1. Upon receipt and review of the Officer's Impaired Driving Affidavit and any associated test results, if applicable, the State of Oklahoma may revoke your Oklahoma driving privileges.
2. If you hold a commercial driver license or were driving a commercial motor vehicle, the State of Oklahoma may disqualify your privilege to operate a commercial motor vehicle.
3. You will continue to have driving privileges, if otherwise eligible, unless or until an Order of Revocation and/or Disqualification is issued by the State of Oklahoma.
4. The State of Oklahoma will send you the Order of Revocation and/or Disqualification to the address you have on file with Service Oklahoma.  
To confirm or change your address on file, contact Service Oklahoma.
5. The Order and Notice of Revocation will inform you of your requirements to regain driving privileges in Oklahoma.

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: \_\_\_\_\_ Place (location when signed): \_\_\_\_\_ Signature of arresting officer: \_\_\_\_\_

|                           |                                |                          |                     |
|---------------------------|--------------------------------|--------------------------|---------------------|
| NAME<br>PAPE, CHRISTOPHER |                                | AGENCY<br>BOARD OF TESTS |                     |
| BADGE<br>0000000484       | AGENCY ADDRESS<br>INTER-AGENCY | ZIP CODE<br>00000        | PHONE<br>4054252460 |



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT  
Chairman


Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

Joshua Smith  
Director

**Certificate of Analysis – Reference Solution**

ALPHA INSTRUMENT  
Intoxilyzer - Alcohol Analyzer  
Model 8000 SN 80-002591  
02/18/2026 11:13

| Test            | g/210L | Time  |
|-----------------|--------|-------|
| Air Blank       | 0.000  | 11:13 |
| Cal Check       | 0.021  | 11:14 |
| Air Blank       | 0.000  | 11:14 |
| Cal Check       | 0.020  | 11:15 |
| Air Blank       | 0.000  | 11:16 |
| Cal Check       | 0.020  | 11:16 |
| Air Blank       | 0.000  | 11:17 |
| Cal Check       | 0.021  | 11:18 |
| Air Blank       | 0.000  | 11:18 |
| Cal Check       | 0.020  | 11:19 |
| Air Blank       | 0.000  | 11:19 |
| Cal Check       | 0.020  | 11:20 |
| Air Blank       | 0.000  | 11:21 |
| Cal Check       | 0.020  | 11:21 |
| Air Blank       | 0.000  | 11:22 |
| Cal Check       | 0.020  | 11:23 |
| Air Blank       | 0.000  | 11:23 |
| Cal Check       | 0.020  | 11:24 |
| Air Blank       | 0.000  | 11:24 |
| Cal Check       | 0.020  | 11:25 |
| Air Blank       | 0.000  | 11:26 |
| Cal Check Stats |        |       |
| Average         | 0.0202 |       |
| Std Dev         | 0.0004 |       |
| Rel Std Dev(%)  | 2.0872 |       |

| NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.   |                                    |                                     |
|--|------------------------------------|-------------------------------------|
| Manufacturer: <b>Guth</b>  | Manufactured Date: <b>9/3/2025</b> | Concentration: <b>0.02</b>          |
| Lot No.: <b>25310</b>  | Expiration Date: <b>9/3/2027</b>   | Equipment Model: <b>Guth 12V500</b> |
| The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration. |                                    |                                     |
| Signature<br>   | Date:<br><b>2/18/2026</b>          |                                     |
| Name <b>Christopher Pape</b>   | Permit No.: <b>23020</b>           |                                     |

**GUTH LABORATORIES, INC.**

690 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-664-6470

**CERTIFICATE OF ANALYSIS**RCVD  
1/20/26  
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **25310** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **September 4, 2025**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0239%** (w/vol) ethyl alcohol. The expiration date for this lot number is **September 3, 2027** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.020 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

***NIST Traceability:***

*Testing was conducted using Cerilliant Reference Standard lot number FN03062417 whose values are traceable to NIST.*

*All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.*



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT  
Chairman

Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

Joshua Smith  
Director

**Certificate of Analysis – Reference Solution**

**ALPHA INSTRUMENT**

Intoxilyzer - Alcohol Analyzer


Model 8000

SN 80-002591

03/30/2026

11:13

| Test            | g/210L | Time  |
|-----------------|--------|-------|
| Air Blank       | 0.000  | 11:14 |
| Cal Check       | 0.040  | 11:14 |
| Air Blank       | 0.000  | 11:15 |
| Cal Check       | 0.039  | 11:16 |
| Air Blank       | 0.000  | 11:16 |
| Cal Check       | 0.039  | 11:17 |
| Air Blank       | 0.000  | 11:18 |
| Cal Check       | 0.040  | 11:18 |
| Air Blank       | 0.000  | 11:19 |
| Cal Check       | 0.039  | 11:20 |
| Air Blank       | 0.000  | 11:20 |
| Cal Check       | 0.040  | 11:21 |
| Air Blank       | 0.000  | 11:21 |
| Cal Check       | 0.039  | 11:22 |
| Air Blank       | 0.000  | 11:23 |
| Cal Check       | 0.040  | 11:23 |
| Air Blank       | 0.000  | 11:24 |
| Cal Check       | 0.039  | 11:25 |
| Air Blank       | 0.000  | 11:25 |
| Cal Check       | 0.039  | 11:26 |
| Air Blank       | 0.000  | 11:26 |
| Cal Check Stats |        |       |
| Average         | 0.0394 |       |
| Std Dev         | 0.0005 |       |
| Rel Std Dev (%) | 1.3106 |       |

| NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.   |                                     |                                     |
|--|-------------------------------------|-------------------------------------|
| Manufacturer: <b>Guth</b>  | Manufactured Date: <b>12/3/2025</b> | Concentration: <b>0.040</b>         |
| Lot No.: <b>25400</b>  | Expiration Date: <b>12/3/2027</b>   | Equipment Model: <b>Guth 12V500</b> |
| The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration. |                                     |                                     |
| Signature<br>   | Date:<br><b>3/30/2026</b>           |                                     |
| Name <b>Christopher Pape</b>   | Permit No.: <b>23020</b>            |                                     |

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

**CERTIFICATE OF ANALYSIS**RCVD  
1/20/26  
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **25400** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **December 4, 2025**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0485%** (w/vol) ethyl alcohol. The expiration date for this lot number is **December 3, 2027** at 11:59 PM.

When used in a calibrated Simulator, operating at  $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$ , this solution will give a breath alcohol analysis instrument reading of **0.040 g/210L**  $\pm 3\%$ .

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

***NIST Traceability:***

*Testing was conducted using Cerilliant Reference Standard lot number FN06032102 whose values are traceable to NIST.*

*All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.*



**STATE OF OKLAHOMA  
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT  
Chairman

Post Office Box 36307  
Oklahoma City, Oklahoma 73136-2307  
Phone: (405) 425-2460 Fax: (405) 425-2490  
www.bot.ok.gov

Joshua Smith  
Director

**Certificate of Analysis – Reference Solution**

ALPHA INSTRUMENT

Intoxilyzer - Alcohol Analyzer


Model 8000

SN 80-002591

03/06/2026

09:29

| Test            | g/210L | Time  |
|-----------------|--------|-------|
| Air Blank       | 0.000  | 09:29 |
| Cal Check       | 0.099  | 09:30 |
| Air Blank       | 0.000  | 09:31 |
| Cal Check       | 0.098  | 09:31 |
| Air Blank       | 0.000  | 09:32 |
| Cal Check       | 0.099  | 09:33 |
| Air Blank       | 0.000  | 09:33 |
| Cal Check       | 0.099  | 09:34 |
| Air Blank       | 0.000  | 09:35 |
| Cal Check       | 0.098  | 09:35 |
| Air Blank       | 0.000  | 09:36 |
| Cal Check       | 0.098  | 09:37 |
| Air Blank       | 0.000  | 09:37 |
| Cal Check       | 0.098  | 09:38 |
| Air Blank       | 0.000  | 09:38 |
| Cal Check       | 0.098  | 09:39 |
| Air Blank       | 0.000  | 09:40 |
| Cal Check       | 0.098  | 09:40 |
| Air Blank       | 0.000  | 09:41 |
| Cal Check       | 0.099  | 09:42 |
| Air Blank       | 0.000  | 09:42 |
| Cal Check Stats |        |       |
| Average         | 0.0984 |       |
| Std Dev         | 0.0005 |       |
| Rel Std Dev (%) | 0.5247 |       |
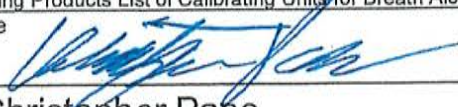
| NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.   |                                      |                                     |
|--|--------------------------------------|-------------------------------------|
| Manufacturer: <b>Guth</b>  | Manufactured Date: <b>08/12/2025</b> | Concentration: <b>0.100</b>         |
| Lot No.: <b>25282</b>  | Expiration Date: <b>08/12/2027</b>   | Equipment Model: <b>Guth 12V500</b> |
| The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration. |                                      |                                     |
| Signature<br>   | Date:<br><b>03/06/2026</b>           |                                     |
| Name<br><b>Christopher Pape</b>  | Permit No.: <b>23020</b>             |                                     |

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

**CERTIFICATE OF ANALYSIS**RCVD  
1/20/26  
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **25282** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **August 14, 2025**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.1209%** (w/vol) ethyl alcohol. The expiration date for this lot number is **August 12, 2027** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

***NIST Traceability:***

*Testing was conducted using Cerilliant Reference Standard lot number FN04042408 whose values are traceable to NIST.*

*All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.*



**STATE OF OKLAHOMA**  
**BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT  
 Chairman

Post Office Box 36307  
 Oklahoma City, Oklahoma 73136-2307  
 Phone: (405) 425-2460 Fax: (405) 425-2490  
 www.bot.ok.gov

Joshua Smith  
 Director

**Certificate of Analysis – Reference Solution**

ALPHA INSTRUMENT

Intoxilyzer - Alcohol Analyzer

Model 8000

SN 80-002591

02/26/2026


09:24

| Test            | g/210L | Time  |
|-----------------|--------|-------|
| Air Blank       | 0.000  | 09:24 |
| Cal Check       | 0.198  | 09:25 |
| Air Blank       | 0.000  | 09:25 |
| Cal Check       | 0.199  | 09:26 |
| Air Blank       | 0.000  | 09:27 |
| Cal Check       | 0.198  | 09:27 |
| Air Blank       | 0.000  | 09:28 |
| Cal Check       | 0.199  | 09:29 |
| Air Blank       | 0.000  | 09:29 |
| Cal Check       | 0.199  | 09:30 |
| Air Blank       | 0.000  | 09:31 |
| Cal Check       | 0.198  | 09:31 |
| Air Blank       | 0.000  | 09:32 |
| Cal Check       | 0.199  | 09:33 |
| Air Blank       | 0.000  | 09:33 |
| Cal Check       | 0.198  | 09:34 |
| Air Blank       | 0.000  | 09:34 |
| Cal Check       | 0.199  | 09:35 |
| Air Blank       | 0.000  | 09:36 |
| Cal Check       | 0.198  | 09:36 |
| Air Blank       | 0.000  | 09:37 |
| Cal Check Stats |        |       |
| Average         | 0.1985 |       |
| Std Dev         | 0.0005 |       |
| Rel Std Dev(%)  | 0.2655 |       |

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.

|                           |                                      |                                     |
|---------------------------|--------------------------------------|-------------------------------------|
| Manufacturer: <b>Guth</b> | Manufactured Date: <b>10/21/2025</b> | Concentration: <b>0.200</b>         |
| Lot No.: <b>25360</b>     | Expiration Date: <b>10/21/2027</b>   | Equipment Model: <b>Guth 12V500</b> |

The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.

|  |                           |
|--|---------------------------|
| Signature<br> | Date:<br><b>2/26/2026</b> |
| Name<br><b>Christopher Pape</b>  | Permit No.: <b>23020</b>  |

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

**CERTIFICATE OF ANALYSIS**RCVD  
1/20/26  
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **25360** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **October 22, 2025**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.2410%** (w/vol) ethyl alcohol. The expiration date for this lot number is **October 21, 2027** at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of **0.200 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

**NIST Traceability:**

Testing was conducted using Certified Reference Standard lot number **FN07012401** whose values are traceable to NIST.  
All balances are calibrated annually by an outside agency using NIST traceable weights.  
Calibration verification is done prior to each use utilizing NIST traceable weights.

3/31/26

Unit Serial #  
80-003314

19/19

Location Returned From  
McCloud PD

Old Cylinder #: 01524080A1-056

Inventory Process

- Unit checked in to BOT
- New Unit Checked out

- Unlink Gas Canister from Unit
- Remove old cylinder from unit (After Lab Evaluation)
- Old cylinder disposed in inventory
- Old cylinder archived

Check in Procedure

- Connect unit to ethernet connection 172.019.004. 137
- Download any stored records 6
- Change unit location to Serial number
- Update Unit Date and time

Previous time: 4/21/26 1101  
 New Time: 4/21/26 1059

Notes from the Field

- Routine Swap
- Issues in the Field

Lab Evaluation

- Breath Hose Hold Vacuum?
- Test Chamber Hold Vacuum?

DVM

3µm 12055 Chamber Temp 47.00  
 9µm 13915 Breath Hose Temp 45

- Snap Test ok?
- RFI Detected?

Peripheral Tests

- Speaker (F5)
- Gas Regulator (F11)
- Power LED Green (F6)
- Pump (F10)
- Red (F7)
- Receipt paper
- Fan (F9)

Atmospheric Sensor

Pressure 973 Correction Factor 1.04

Flow Sensor Resistance 236

Tank Sensor Value 939

Is regulator closed?  Yes /  No

Sent to CMI (If needed)

Date Sent \_\_\_\_\_  
 Invoice Sent \_\_\_\_\_  
 Date Returned \_\_\_\_\_

Reason Sent

Repairs needed:

Repairs Made:

Ready for Bench Check?

- To complete if unit sent to CMI
- CMI Workorder
- CMI Certificate of Calibration
- CMI Service Evaluation Form

Inventory Process

- Linked Cylinder to Unit
- Changed Unit Status
- Lab Check Report Completed

- Bench Check Technician Report
- Mock Subject Test Completed
- Certificate of Calibration and Operation
- Service Overview Completed

For BOT Internal Use Only

Intoxily