



**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Post Office Box 36307

Oklahoma City, Oklahoma 73136-2307

Phone: (405) 425-2460 Fax: (405) 425-2490

www.bot.ok.gov

Dr. Jarrad Wagner, Ph.D. F-ABFT
Chairman

Joshua Smith
Director

ATTESTATION

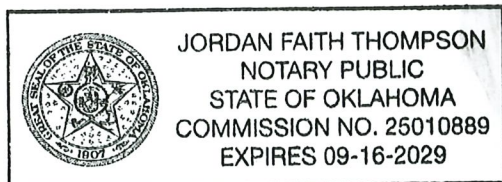
I, the undersigned Joshua Smith (Director/Records Custodian), attest under oath pursuant to 12 O.S. §2902 this record was made, at or near the time of the occurrence of the matters set forth by or from information transmitted by a person having knowledge of those matters; was kept in the course of the regularly conducted business activity of the Board; and was made pursuant to the regularly conducted activity. The duplicate record provided is listed below:

- 1) The attached "Service Overview" is a true and accurate copy of maintenance performed on the Intoxilyzer 8000, serial number 80-006160, in compliance with the agency's Maintenance Policy and Procedures.

Signed _____

Joshua Smith

Name of Position: Director/Records Custodian



State of Oklahoma

County of Oklahoma

Signed or attested before me on 10/01/2025 by

Joshua Smith.

Jordan Faith Thompson
(Signature of notarial officer)

Service Overview

Board of Tests for Alcohol and Drug Influence
*Enhancing public safety by ensuring the accuracy and scientific reliability
of blood and breath testing.*



2/20

3600 N. Martin Luther King Bldg #9
Oklahoma City, OK 73111
405-425-2460

BENCH CHECK DATE: 06/13/2022

SERIAL: 80-006160

DOCUMENTATION WITHIN SERVICE PACKET:

- BOT Certificate of Maintenance, Calibration and Operation
- Bench Check
- CMI Certificate of Calibration
- CMI Evaluation Form
- CMI Work Order
- BOT Dry Gas Certificate of Analysis
- Ilmo Specialty Gases Certificate of Analysis
- Mock Subject Test
- BOT Certificate of Analysis - Reference Solution (4)
- Guth Laboratories Certificate of Analysis (4)



3/20

STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

CERTIFICATE of MAINTENANCE, CALIBRATION, and OPERATION


This certifies the **Intoxilyzer 8000** referenced below, manufactured by CMI, Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, is maintained by the Board of Tests for Alcohol and Drug Influence in accordance with 40 O.A.C. 30-1-3. The instrument reported below is hereby deemed to be in compliance with Board rules via linearity accuracy examination and operational/functional procedure evaluation also referred to as a "Bench Check". The Intoxilyzer 8000 appears on the National Highway Traffic Safety Administration's Conforming Products List of Evidential Breath Alcohol Measurement Devices (82 FR 50940 *et seq*) as an approved device. In accordance with 40 O.A.C. 25-1-2, this instrument is approved and certified for use in the State of Oklahoma as an automated evidential breath-alcohol measurement device by the Board of Tests pursuant to 47 O.S. § 759.

All analyses conducted by this instrument are considered precise, valid, and in accordance with 40 O.A.C. 30-1-3(b) Operating Procedure(s). A printed breath test operator affidavit containing the seal of the Board by this instrument confirms a calibration check successfully verified the instrument's calibration at the time of test. The printed breath test affidavit also serves as proof that safeguards in place to protect the breath samples, including but not limited to Ambient Air Interference, Radio Frequency Interference, Residual Mouth Alcohol, Deficient Samples, etc. were not triggered and the test is a true and accurate scientific measurement of the person's alcohol concentration.

- Maintenance was performed on the Intoxilyzer 8000, Serial No. 80-006160 on 06/13/2022 in accordance with 40 O.A.C. 30-1-3 and internal BOT maintenance policy and procedures.
- A dry gas canister appearing on the National Highway Traffic Safety Administration's Conforming Products List of Calibrating Units for Breath Alcohol Testers (77 FR 64588 *et seq*) and approved in accordance with 40 O.A.C. 25-1-3 and BOT maintenance policy was installed.
 - Manufacturer: ILMO Unit: 208 ppm / 105 L
 - Cylinder Lot No. - Canister No.: 06022080A2-051
 - BOT analysis average test result: 0.081
 - Cylinder expiration date: 04/05/2024
- In accordance with 40 O.A.C. 25-1-2.1, the approved reference standards and equipment were used:

Simulator Model			
Guth 2100	Guth 2100	Guth 2100	Guth 2100
Serial #			
DR3685	DR3694	DR3594	DR3754
Concentration			
.02	.04	.10	.20
LOT #			
21120	21070	21190	20510
Manufactured Date			
4/7/2021	3/1/2021	6/8/2021	12/3/2020
Expiration Date			
4/7/2023	3/1/2023	6/8/2023	12/3/2022
Solution Commission Date			
06/08/2022	06/08/2022	06/08/2022	06/08/2022

Completed by BOT Personnel:


Signature/Name Christopher Pape

Breath-Alcohol Testing Program Administrator
Title


Signature Name Joshua Smith

State Director
Title

STATE OF OKLAHOMA - BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE

CMI INTOXILYZER BENCH CHECK PROCEDURE REPORT

BENCH CHECK DATE 06/13/2022	TIME PROCEDURE BEGAN 11:03:4	INSTRUMENT LOCATION 80-006160
INSTRUMENT TYPE Intoxilyzer	INSTRUMENT MODEL I-8000	INSTRUMENT SERIAL NUMBER 80-006160

TEST RESULTS

Diagnostics**PASS****Mouth Alcohol Test****PASS****RFI Test****PASS****Abort Test****PASS****Improper Sample Test****PASS****Interferent Detect Test****PASS****No Sample Given Test****PASS****Insufficient Sample Test****PASS****Air Blank Contamination Test****PASS****Range Exceeded Test****PASS**

Duplicate Breath Sampling Results

Value	Test 1	Test 2	Control Test
.02	0.019	0.019	0.020
.04	0.039	0.038	0.040
.10	0.100	0.099	0.101
.20	0.201	0.198	0.197

PASS

CERTIFYING TECHNICIAN

I swear under penalty of perjury that in regards to the above listed instrument, I have complied with the Intoxilyzer Bench Checking Procedure approved by the State of Oklahoma.

SIGNATURE 	DATE 6/13/2022
NAME PAPE, CHRISTOPHER	PERMIT # 0000023020

Certificate of Calibration

This is to certify the calibration of **Intoxilyzer**® serial number 80-006160, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58, No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date June 2, 2022

Signed

[Signature]

Technician

EMI INC.

316 East 9th Street
Owensboro, KY 42303
USA

Part No. 650517 Rev.A

CMI INC.

Service Evaluation Form

This form **MUST** be completed and enclosed with item to be serviced.
 *** **Hazardous Material Warning! – DO NOT return gas cylinder with instrument!** ***

Note: Please ship items in their original shipping container.

Contact information:

Name Joshua Smith **Phone:** (405) 425-2460
Email: joshua.smith@bot.ok.gov **Customer #** 731362 (contact Customer Service)

Your Billing Address

Agency Name: Oklahoma Board of Tests for Alcohol and Drug Influence

Address: PO Box 36307

City: Oklahoma City **State:** OK **Zip:** 73136

Credit Card/PO #:

Name on Card:

Expiration Date: _____ **CVV** _____

Your Shipping Address

Agency Name: Oklahoma Board of Tests for Alcohol and Drug Influence

Address: 3600 N. Martin Luther King Ave

City: Oklahoma City **State:** OK **Zip:** 73111

Contact Person: Christopher Pape

Contact Phone: 405-425-2256

Email: Christopher.Pape@bot.ok.gov

Instrument Serial Number: 80-006160

Detailed Description of Problem:

Unit will not post after being powered. (NO BREATH HOSE) May have records in memory.

Note: For instruments not under warranty, an evaluation fee of \$83.⁰⁰ (infrared and I-200) or \$43.⁵⁰ (fuel cell), plus return shipping cost, will apply to all service items.

☒ **I authorize all repairs to be performed.**
 Any repairs not exceeding \$1,500 require no authorization.

After repair, add "Certificate of Calibration" for: ☐ S-D2, S-D5, and I-800 (\$19)

☐ I-200, I-240, I-300, I-400, and I-600 (\$39)

Authorized By:

Joshua Smith
 Name (Please Print)

Director

Title

Signature

03/01/2022
 Date

Ship item to:

CMI, Inc.
Attn: Service Dept.
316 East Ninth Street
Owensboro, KY 42303



7/20

DATE: 3/9/2022

316 E 9th Street / Owensboro KY 42303 / USA
Phone: 866-835-0690 Fax: 270-685-6268

SERVICE WORK ORDER

406765

Bill To:

Oklahoma Board Of Tests
Blvd Adcu Annex Bldg #9
3600 Martin Luther King
Oklahoma City, OK 73136-3000
USA
Customer #: 731362
Phone: 405-425-2424
Fax: 405-425-2490 fax

Ship To:

Oklahoma Board Of Tests
Blvd Adcu Annex Bldg #9
3600 Martin Luther King
Oklahoma City OK 73136-3000
USA
Phone: 405-425-2424

MODEL #: 002480OK

Serial Number: 80-006160

BILL CODE: Out of Warranty

EXTRA PARTS RCVD:

KB,BATTERY COVER,GAS SHELF W.REG (no bh)

DESCRIPTION OF PROBLEM

UNIT WILL NOT POST AFTER BEING POWERED. (NO BREATH HOSE) MAY HAVE RECORDS IN MEMORY.

WORK PERFORMED:

No Breath Hose; Data Flash/Corrup; Prom Flash/Corrupt; Replaced all Parts listed; Calibrated and Final Tested;

PARTS USED

Seq. No.	Part	Description	Quantity
100	650517	CERTIFICATE OF CALIBRATION	2.00 EA
110	441169	COVER DUST,5/8in x 1/2in	1.00 EA
120	470145	CAP,PLASTIC,.25IDx.50LG,RED	3.00 EA
130	470154	CAP,PLSTC,.406IDx.438-.562,RED	1.00 EA
135	220284	IC,SMD,FLASH,512KB,PLCC-32	1.00 EA
136	021307FL	HOSE ASSY,BREATH,I8000	1.00 EA
137	230115	IC,SMD,DATAFLASH,4MB,TSOP-28	2.00 EA
140	470218	SEAL,TAMPER EVIDENT,SERVICE	1.00 EA

LABOR / TESTING

Misc Code	Description	Hours
LABR	Service Repair Labor	
LCAL	Service - Calibration Adjust	0.50
LFT	Service - Cal / Final Test	1.00

Technician Name: Hayden Schweikardt

SRV4067650001





8/20

**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

Certificate of Analysis – Dry Gas

ALPHA INSTRUMENT

Intoxilyzer – Alcohol Analyzer

Model 8000

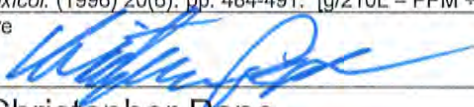
SN 80-002591

04/25/2022

12:39

Test	g/210L	Time
Air Blank	0.000	12:39
Cal Check	0.081	12:40
Air Blank	0.000	12:40
Cal Check	0.081	12:40
Air Blank	0.000	12:41
Cal Check	0.080	12:41
Air Blank	0.000	12:42
Cal Check	0.081	12:42
Air Blank	0.000	12:43
Cal Check	0.081	12:43
Air Blank	0.000	12:44
Cal Check	0.081	12:44
Air Blank	0.000	12:44
Cal Check	0.081	12:45
Air Blank	0.000	12:45
Cal Check	0.080	12:46
Air Blank	0.000	12:46
Cal Check	0.080	12:46
Air Blank	0.000	12:47
Cal Check	0.081	12:47
Air Blank	0.000	12:48
Cal Check Stats		
Average	0.0807	
Std Dev	0.0004	
Rel Std Dev (%)	0.5985	

Known Breath-Alcohol Equivalent (g/210L) Ratio of Nitrogen and Ethanol Gas (208 ppm/105 L)
0.08 g/210L = 208 PPM ÷ 2605*

Manufacturer: ILMO	Manufacturer Cert. No.: 14689	Cylinder Size: 105L
Lot No.: 06022080A2	Cylinder No.: 51	Expiration Date: 04/05/2024
The above 105 L cylinder has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference standard meets the requirements for use as a calibrating unit in accordance with 47 O.S. 759, 40 OAC 25-1-3, and appears on the current conforming products list of calibrating units for breath alcohol testers published by the U.S. Department of Transportation in the Federal Register [77 FR 64588 et seq]. * See Dubowski, K.M. and Essary N., <i>Journal of Anal. Toxicol.</i> (1996) 20(6): pp. 484-491. [g/210L = PPM ÷ 2605].		
Signature: 		Date: 04/25/2022
Name: Christopher Pape		Permit No.: 23020



7 Eastgate Dr. • P.O. Box 790 • Jacksonville, IL 62651-0790
217-245-2183 • Fax: 217-243-7634 • www.ilmoproducts.com

Certificate of Analysis

Certificate ID: 14689 ✓
Part #: BAC105L080T
Cylinder Size: 105L
Lot Number: 06022080A2 ✓
Expiration: 4/5/2024 ✓

RCV'D
4/11/2022
CP

0.080 BAC (For the calibration of instruments used to determine breath alcohol concentration)

Contents: 105 Liters @ 1000 psig 70°F (21°C)

Component:	Reported Concentration:	Analytical Accuracy (U, k=2):	Analytical Method:
Ethanol	208 ppm	+/-0.002 BAC(6/210L)	NDIR
Nitrogen	Balance	[5.2 ppm]	

Distributed by:

CMI Inc.
316 East Ninth Street
Owensboro, KY 42303
Phone 866-835-0690
www.alcoholtest.com

*Traceable to:

Certified Reference Material - 251.8 µmol/mol

Ethanol in Nitrogen - Serial No. CC198518 Lot No. 111318E1

Store in dry area, away from sources of heat, ignition and direct sunlight. Do not allow storage area to exceed 52 °C (125 °F).

Specialty Gas Lab Tech

03-28-2022
Issuance Date



The calibration results within this certificate were obtained using equipment and standards capable of producing analytical results traceable to NIST, and apply only to the items contained on this certificate. ILMO Products Company makes no warranty or representation as to the suitability of the use of any information provided for any particular purpose. The information use is at the sole discretion and risk of the user. Liability shall be limited to established replacement cost of this material or service.

ISO/IEC 17025:2017 Accredited Laboratory

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

10/20

Section 1

ARREST DATE 06/13/2022	ARREST TIME 1200	COUNTY # 55	CITY #	CITATION #
ARREST LOCATION BOT		CITY	COUNTY OKLAHOMA	
SUBJECT NAME BREATHTEST, DELTA		DATE OF BIRTH 01/01/1985	SEX M	HEIGHT 510
ADDRESS 3600 MLK		CITY OKLAHOMA CITY	STATE OK	WEIGHT 200
DRIVER LICENSE # B083463042		EXPIRATION DATE 10/31/2022	CLASS D	RESTRICTIONS 73111
VEHICLE MAKE MOCK		MODEL MOCK	TAG # MOCK	STATE OK
			CMV N	PLACARD N

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST


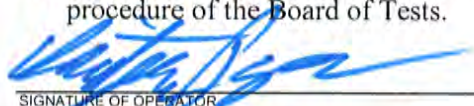
Describe person's condition (odor, actions, etc.):

POST REPAIR MAINTENANCE, INSTALLED NEW CYLINDER

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

EQUIPMENT I-8000	SERIAL NUMBER 80-006160	SPECIALIST BOARD OF TESTS																																			
STANDARD Dry	TARGET VALUE 0.080	PERMIT NUMBER 23020																																			
MANUFACTURER ILMO	LOT # 06022080A2	DATE CYLINDER INSTALLED 06/13/2022	DATE CYLINDER EXPIRES 04/05/2024																																		
OBSERVATION START 1210	TEST DATE 06/13/2022	A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.																																			
RESULT TABLE																																					
<table><thead><tr><th>Test</th><th>g/210L</th><th>Time</th></tr></thead><tbody><tr><td>Air Blank</td><td>0.00</td><td>12:37</td></tr><tr><td>Subject Test</td><td>0.00</td><td>12:37</td></tr><tr><td>Breath Volume</td><td>1.667 LITERS</td><td></td></tr><tr><td>Air Blank</td><td>0.00</td><td>12:38</td></tr><tr><td>Wait</td><td></td><td>12:40</td></tr><tr><td>Air Blank</td><td>0.00</td><td>12:40</td></tr><tr><td>Subject Test</td><td>0.00</td><td>12:41</td></tr><tr><td>Breath Volume</td><td>1.894 LITERS</td><td></td></tr><tr><td>Air Blank</td><td>0.00</td><td>12:41</td></tr><tr><td>Cal Check</td><td>0.08</td><td>12:41</td></tr><tr><td>Air Blank</td><td>0.00</td><td>12:42</td></tr></tbody></table>				Test	g/210L	Time	Air Blank	0.00	12:37	Subject Test	0.00	12:37	Breath Volume	1.667 LITERS		Air Blank	0.00	12:38	Wait		12:40	Air Blank	0.00	12:40	Subject Test	0.00	12:41	Breath Volume	1.894 LITERS		Air Blank	0.00	12:41	Cal Check	0.08	12:41	Air Blank
Test	g/210L	Time																																			
Air Blank	0.00	12:37																																			
Subject Test	0.00	12:37																																			
Breath Volume	1.667 LITERS																																				
Air Blank	0.00	12:38																																			
Wait		12:40																																			
Air Blank	0.00	12:40																																			
Subject Test	0.00	12:41																																			
Breath Volume	1.894 LITERS																																				
Air Blank	0.00	12:41																																			
Cal Check	0.08	12:41																																			
Air Blank	0.00	12:42																																			
		I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.																																			
		SIGNATURE OF OPERATOR 																																			
		NAME PAPE, CHRISTOPHER																																			
		BADGE # 0000000484																																			
		PERMIT # 0000023020																																			
		AGENCY BOARD OF TESTS																																			

Test Result BrAC g/210L 0.00

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. You may appeal any Departmental action against your driver license as a result of this arrest by filing a petition in the District Court of the County in which you were arrested. Your petition must be filed within thirty (30) days of the date listed below.

Section 3

You may be eligible to participate in the Impaired Driver Accountability Program (IDAP) administered by the Department of Public Safety (DPS). Completion of IDAP may prevent a revocation as a result of this arrest from appearing on your driving record. Participation in IDAP may reduce the amount of fees you will be required to pay to the State. You must request IDAP within thirty (30) days of the date listed below. Mail the request for IDAP to the Department of Public Safety at P.O. Box 11415, Oklahoma City, OK 73136.

SERVING OF NOTICE: I personally hand delivered a copy of the above containing the Notice of Revocation and the Notice of Department Impaired Driving Accountability Program to the person named in Section 1 above.

Date served: _____ Signature of officer: _____ Badge #: _____

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: _____ Place (location when signed): _____ Signature of arresting officer: _____

NAME PAPE, CHRISTOPHER	AGENCY BOARD OF TESTS
BADGE 0000000484	AGENCY ADDRESS INTER-AGENCY
ZIP CODE 00000	PHONE 4054252460

OTHER WITNESSES:

Section 5

1	NAME	TITLE	ADDRESS	PHONE
2	NAME	TITLE	ADDRESS	PHONE

OFFICER/DPS COPY

BOT Form AFF01 09/2021

OFFICER'S AFFIDAVIT AND NOTICE OF REVOCATION/DISQUALIFICATION

11/20

Section 1

ARREST DATE 06/13/2022	ARREST TIME 1200	COUNTY # 55	CITY #	CITATION #
ARREST LOCATION BOT		CITY	COUNTY OKLAHOMA	
SUBJECT NAME BREATHTEST, DELTA		DATE OF BIRTH 01/01/1985	SEX M	HEIGHT 510
ADDRESS 3600 MLK		CITY OKLAHOMA CITY	STATE OK	WEIGHT 200
DRIVER LICENSE # B083463042		EXPIRATION DATE 10/31/2022	CLASS D	RESTRICTIONS
VEHICLE MAKE MOCK		MODEL MOCK	TAG # MOCK	ENDORSEMENTS
		STATE OK	CMV N	PLACARD N

On the above date, time, and location, the above named person was arrested, and I had reasonable grounds to believe the person was driving or in actual physical control of a motor vehicle upon a public road, highway, street, turnpike, public parking lot, or other public place, or upon a private road, street, alley or lane which provides access to one or more single or multi-family dwelling within this state while under the influence of alcohol and/or other intoxicating substances as prohibited by law.

(Describe driving behavior or circumstances):

MOCK SUBJECT TEST


Describe person's condition (odor, actions, etc.):

POST REPAIR MAINTENANCE, INSTALLED NEW CYLINDER

THE PERSON WAS INFORMED OF THE IMPLIED CONSENT TEST REQUEST

BREATH-ALCOHOL ANALYSIS RECORD, REPORT OF THE PERSON NAMED IN SECTION 1 ABOVE, and LOG OF TEST AND MAINTENANCE RECORD.

Section 2

EQUIPMENT I-8000	SERIAL NUMBER 80-006160	SPECIALIST BOARD OF TESTS	
STANDARD Dry	TARGET VALUE 0.080	PERMIT NUMBER 23020	
MANUFACTURER ILMO	LOT # 06022080A2	DATE CYLINDER INSTALLED 06/13/2022	DATE CYLINDER EXPIRES 04/05/2024
OBSERVATION START 1210	TEST DATE 06/13/2022	A 15-minute (or greater) deprivation/observation period was observed prior to and continuing to the time of test administration.	
RESULT TABLE		I administered said Breath-Alcohol Test in compliance with the applicable operating procedure of the Board of Tests.	
Test	g/210L	Time	
Air Blank	0.00	12:37	
Subject Test	0.00	12:37	
Breath Volume	1.667 LITERS		
Air Blank	0.00	12:38	
Wait		12:40	
Air Blank	0.00	12:40	
Subject Test	0.00	12:41	
Breath Volume	1.894 LITERS		
Air Blank	0.00	12:41	
Cal Check	0.08	12:41	
Air Blank	0.00	12:42	
NAME PAPE, CHRISTOPHER			Test Result BrAC g/210L 0.00
BADGE # 0000000484			
PERMIT # 0000023020			
AGENCY BOARD OF TESTS			

EFFECTIVE THIRTY (30) DAYS FROM DATE SERVED. Notice of Revocation is given that your driver's license is revoked or denied for a period of one hundred eighty (180) days or more, if you refused or failed the state's test(s). In addition, your commercial driving privilege in this state will be disqualified if you refused or failed the state's test(s) (operating a CMV or CDL holder) OR a test result of .04 or more while operating a commercial motor vehicle. You may appeal any Departmental action against your driver license as a result of this arrest by filing a petition in the District Court of the County in which you were arrested. Your petition must be filed within thirty (30) days of the date listed below.

You may be eligible to participate in the Impaired Driver Accountability Program (IDAP) administered by the Department of Public Safety (DPS). Completion of IDAP may prevent a revocation as a result of this arrest from appearing on your driving record. Participation in IDAP may reduce the amount of fees you will be required to pay to the State. You must request IDAP within thirty (30) days of the date listed below. Mail the request for IDAP to the Department of Public Safety at P.O. Box 11415, Oklahoma City, OK 73136.

SERVING OF NOTICE: I personally hand delivered a copy of the above containing the Notice of Revocation and the Notice of Department Impaired Driving Accountability Program to the person named in Section 1 above.

Date served: _____ Signature of officer: _____ Badge #: _____

In accordance with Title 12 O.S. Section 426, "I state under penalty of perjury under the laws of Oklahoma that the foregoing is true and correct."

Section 4

Date: _____	Place (location when signed): _____	Signature of arresting officer: _____
NAME PAPE, CHRISTOPHER	AGENCY BOARD OF TESTS	
BADGE 0000000484	AGENCY ADDRESS INTER-AGENCY	ZIP CODE 00000
		PHONE 4054252460



12/20

**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

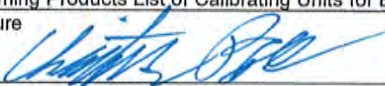
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT 2
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-004923
03/08/2022 11:33

Test	g/210L	Time
Air Blank	0.000	11:34
Cal Check	0.022	11:34
Air Blank	0.000	11:35
Cal Check	0.021	11:36
Air Blank	0.000	11:36
Cal Check	0.021	11:37
Air Blank	0.000	11:37
Cal Check	0.021	11:38
Air Blank	0.000	11:38
Cal Check	0.021	11:39
Air Blank	0.000	11:40
Cal Check	0.021	11:40
Air Blank	0.000	11:41
Cal Check	0.020	11:42
Air Blank	0.000	11:42
Cal Check	0.021	11:43
Air Blank	0.000	11:43
Cal Check	0.021	11:44
Air Blank	0.000	11:45
Cal Check	0.020	11:45
Air Blank	0.000	11:46
Cal Check Stats		
Average	0.0209	
Std Dev	0.0005	
Rel Std Dev(%)	2.7158	

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.		
Manufacturer: Guth	Manufactured Date: 04/07/2021	Concentration: 0.02
Lot No.: 21120	Expiration Date: 04/07/2023	Equipment Model: Guth 2100
The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.		
Signature 	Date: 03/08/2022	
Name Christopher Pape	Permit No.: 23020	



13/20

**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

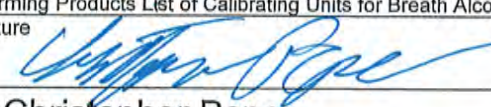
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT 2
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-004923
03/08/2022 11:53

Test	g/210L	Time
Air Blank	0.000	11:53
Cal Check	0.041	11:54
Air Blank	0.000	11:54
Cal Check	0.041	11:55
Air Blank	0.000	11:56
Cal Check	0.042	11:56
Air Blank	0.000	11:57
Cal Check	0.042	11:58
Air Blank	0.000	11:58
Cal Check	0.042	11:59
Air Blank	0.000	11:59
Cal Check	0.041	12:00
Air Blank	0.000	12:01
Cal Check	0.041	12:01
Air Blank	0.000	12:02
Cal Check	0.041	12:02
Air Blank	0.000	12:03
Cal Check	0.041	12:04
Air Blank	0.000	12:04
Cal Check	0.042	12:05
Air Blank	0.000	12:05
Cal Check Stats		
Average	0.0414	
Std Dev	0.0005	
Rel Std Dev(%)	1.2473	

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.		
Manufacturer: Guth	Manufactured Date: 03/01/2021	Concentration: 0.04
Lot No.: 21070	Expiration Date: 03/01/2023	Equipment Model: Guth 2100
The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.		
Signature 	Date: 03/08/2022	
Name Christopher Pape	Permit No.: 23020	



14/20

**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

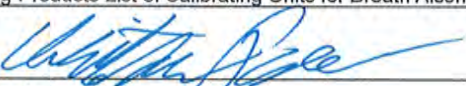
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT 2
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-004923
03/08/2022 12:11

Test	g/210L	Time
Air Blank	0.000	12:11
Cal Check	0.101	12:12
Air Blank	0.000	12:13
Cal Check	0.101	12:13
Air Blank	0.000	12:14
Cal Check	0.101	12:15
Air Blank	0.000	12:15
Cal Check	0.102	12:16
Air Blank	0.000	12:16
Cal Check	0.101	12:17
Air Blank	0.000	12:18
Cal Check	0.101	12:18
Air Blank	0.000	12:19
Cal Check	0.101	12:20
Air Blank	0.000	12:20
Cal Check	0.101	12:21
Air Blank	0.000	12:21
Cal Check	0.101	12:22
Air Blank	0.000	12:23
Cal Check	0.101	12:23
Air Blank	0.000	12:24
Cal Check Stats		
Average	0.1011	
Std Dev	0.0003	
Rel Std Dev(%)	0.3127	

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.		
Manufacturer: Guth	Manufactured Date: 06/08/2021	Concentration: 0.10
Lot No.: 21190	Expiration Date: 06/08/2023	Equipment Model: Guth 2100
The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.		
Signature 	Date: 03/08/2022	
Name Christopher Pape	Permit No.: 23020	



15/20

**STATE OF OKLAHOMA
BOARD OF TESTS FOR ALCOHOL AND DRUG INFLUENCE**

Dr. Jarrad Wagner, Ph.D., F-ABFT
Chairman

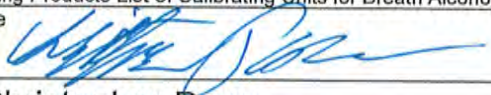
Post Office Box 36307
Oklahoma City, Oklahoma 73136-2307
Phone: (405) 425-2460 Fax: (405) 425-2490
www.bot.ok.gov

Joshua Smith
Director

Certificate of Analysis – Reference Solution

ALPHA INSTRUMENT 2
Intoxilyzer - Alcohol Analyzer
Model 8000 SN 80-004923
03/08/2022 12:28

Test	g/210L	Time
Air Blank	0.000	12:29
Cal Check	0.201	12:29
Air Blank	0.000	12:30
Cal Check	0.200	12:31
Air Blank	0.000	12:31
Cal Check	0.201	12:32
Air Blank	0.000	12:33
Cal Check	0.201	12:33
Air Blank	0.000	12:34
Cal Check	0.201	12:35
Air Blank	0.000	12:35
Cal Check	0.201	12:36
Air Blank	0.000	12:36
Cal Check	0.201	12:37
Air Blank	0.000	12:38
Cal Check	0.200	12:38
Air Blank	0.000	12:39
Cal Check	0.201	12:40
Air Blank	0.000	12:40
Cal Check	0.200	12:41
Air Blank	0.000	12:41
Cal Check Stats		
Average	0.2007	
Std Dev	0.0004	
Rel Std Dev(%)	0.2406	

NIST Traceable Breath-Alcohol Equivalent (g/210L +/- 3%) Ratio of ethanol and water free of test interfering substances.		
Manufacturer: Guth	Manufactured Date: 12/03/2020	Concentration: 0.20
Lot No.: 20510	Expiration Date: 12/03/2022	Equipment Model: Guth 2100
The above solution has been analyzed via Infrared Spectroscopy utilizing the Intoxilyzer 8000 and was performed in accordance with the procedures prescribed by the State Director of Tests pursuant to 40 OAC 30-1-3. The above NIST traceable reference solution meets the requirements for use pursuant to 40 OAC 25-1-2.1. Approved breath-alcohol measurement equipment and reference solutions. Equipment used appears on the current or supplemented Conforming Products List of Calibrating Units for Breath Alcohol Testers (72 FR 34747), published by the National Highway Traffic Safety Administration.		
Signature 	Date: 03/08/2022	
Name: Christopher Pape	Permit No.: 23020	

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **21120** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **April 16, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0237%** (w/vol) ethyl alcohol. The expiration date for this lot number is **April 7, 2023** at 11:59 PM.

When used in a calibrated Simulator, operating at **34°C +/- .2°C**, this solution will give a breath alcohol analysis instrument reading of **0.020 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

*Testing was conducted using Cerilliant Reference Standard lot number **FN08031603** whose values are traceable to NIST.*

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

RCV'D
8/26/21
CP

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSISRCV'D
8/26/21
CP

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **21070** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **March 2, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.0488%** (w/vol) ethyl alcohol. The expiration date for this lot number is **March 1, 2023** at 11:59 PM.

When used in a calibrated Simulator, operating at $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$, this solution will give a breath alcohol analysis instrument reading of **0.040 g/210L \pm 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

*Testing was conducted using Cerilliant Reference Standard lot number **FN08211802** whose values are traceable to NIST.*

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSISRCV'D
8/26/21
CP**Certified Alcohol Reference Solution for Simulator**

Random Samples of Lot Number **21190** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on **June 9, 2021**, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain **0.1206%** (w/vol) ethyl alcohol. The expiration date for this lot number is **June 8, 2023** at 11:59 PM.

When used in a calibrated Simulator, operating at **34°C +/- .2°C**, this solution will give a breath alcohol analysis instrument reading of **0.100 g/210L +/- 3%**.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN03052002 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number 20510 of Alcohol Reference Solution for Simulator were analyzed by gas chromatography on December 7, 2020, using a Perkin Elmer Gas Chromatograph Autosystem XL S/N: 610N9030209, and found to contain 0.2419% (w/vol) ethyl alcohol. The expiration date for this lot number is December 3, 2022 at 11:59 PM.

When used in a calibrated Simulator, operating at 34°C +/- .2°C, this solution will give a breath alcohol analysis instrument reading of 0.200 g/210L +/- 3%.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President
GUTH LABORATORIES, INC.

NIST Traceability:

Testing was conducted using Cerilliant Reference Standard lot number FN06231703 whose values are traceable to NIST.

All balances are calibrated annually by an outside agency using NIST traceable weights. Calibration verification is done prior to each use utilizing NIST traceable weights.

RCV'D
8/26/21
CP

Location Returned From Owasso PD	Unit Serial # 80-006160	Date: 20/20 2/1/22
-------------------------------------	----------------------------	-----------------------

Inventory Process

☒ Unit checked in to BOT
 ☒ New Unit Checked out
 ☒ Unlink Gas Canister from Unit
 ☐ Remove old cylinder from unit (After Lab Evaluation)
 ☐ Old cylinder archived
 ☐ Old cylinder disposed in inventory
 ☐ Old cylinder archived

Old Cylinder #: 22620080A2-002

<p><u>Check in Procedure</u></p> <p> <input type="checkbox"/> Connect unit to ethernet connection 172.019.004. 133 <input type="checkbox"/> Download any stored records _____ <input type="checkbox"/> Change unit location to Serial number <input type="checkbox"/> Update Unit Date and time Previous time: _____ New Time: _____ <input type="checkbox"/> Execute a new Forms load on the unit </p>	<p><u>Notes from the Field</u></p> <p> <input type="checkbox"/> Routine Swap <input checked="" type="checkbox"/> Issues in the Field </p> <p>Unit will not POST</p>
---	---

Lab Evaluation

Breath Hose Hold Vacuum?: ☐ Test Chamber Hold Vacuum?: ☐

DVM

3µm _____ Chamber Temp _____ Snap Test ok? ☐
 9µm _____ Breath Hose Temp _____ RFI Detected? ☐

Peripheral Tests

Speaker (F5) ☐ Power LED Green (F6) ☐ Red (F7) ☐ Fan (F9) ☐
 Gas Regulator (F11) ☐ Pump (F10) ☐ Receipt paper ☐

Atmospheric Sensor

Pressure _____ Correction Factor _____ Resistance _____ Value _____
 Flow Sensor _____ Tank Sensor _____
 Is regulator close to Tank Sensor Value? Yes / No

Sent to CMI (If needed) Reason Sent

Date Sent 03/01/2022 unit will not post. no breath hose.

Invoice Sent _____

Date Returned _____

Repairs needed:	Repairs Made:
	Ready for Bench Check? <input type="checkbox"/>

<p><u>Inventory Process</u></p> <p> <input type="checkbox"/> Linked Cylinder to Unit <input type="checkbox"/> Changed Unit Status <input type="checkbox"/> Bench Check Report Completed </p>	<p> <input type="checkbox"/> Bench Check Technician Report <input type="checkbox"/> Mock Subject Test Completed <input type="checkbox"/> Certificate of Calibration and Operation <input type="checkbox"/> Service Overview Completed </p>	<p>To complete if unit sent to CMI</p> <p> <input type="checkbox"/> CMI Workorder <input type="checkbox"/> CMI Certificate of Calibration _____ <input type="checkbox"/> CMI Service Evaluation Form </p>
--	---	---